



Lympha Press USA

LYMPHAPOD®

Features:

- Lympha Press® calibrated gradient compression
- Treats legs, abdomen, hips and genital areas
- 40 overlapping air chambers
- Fully adjustable pressure in 4 - 12 zones
- Designed for large patients with limited mobility
- One size fits most patients from 350 to 700 lbs.



A sequential pneumatic compression therapy system designed for the larger patient

Larger patients with lymphedema have special needs. They find it difficult to manage bandaging, self massage and other treatment systems because of limited mobility and the large size of the affected limb. They may also have severe lymphedema on the upper thighs and torso, which cannot be treated with standard leg garments.

LymphaPod® is a home care system designed for ease of use by large patients with limited mobility. It treats the entire lower body, including legs, hips, and abdomen.

Ways to use LYMPHAPOD®:

At home: Pneumatic compression works best with daily use. Patients can benefit most from treatments at home. It is also cost effective.

In clinic: LymphaPod® is used in hospitals and clinics as part of a program of complete decongestive therapy.

Use LYMPHAPOD® with:




Lympha Press Optimal®

Pressure and treatment time are all adjustable to meet individual needs, with patient lock and easy on/off operation for home use. Medicare coded E0652.

The LymphaPod® treatment system is designed specifically for larger patients with lymphedema or venous insufficiency

- One size fits most patients from 350-700+ lbs. with ease.
- Treats entire lower body.
- Special emphasis on the inner thigh, to treat lobules.
- Ergonomically designed for patients with limited mobility.*
- Lympha Press® sequential massage, with fully adjustable pressure in 4-12 zones.


* Patients with limited mobility may need the assistance of a caregiver.

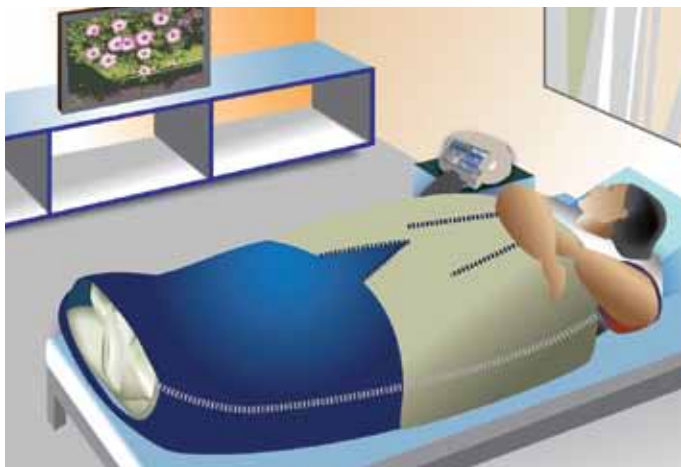

LymphaPod® dimensions in cm.

LymphaPod® Left Hose Exit Cat. No. L4MP800L		LymphaPod® Right Hose Exit Cat. No. L4MP800R		
Reference Point	On patient (approximate*)	Circumference (cm)	With one expander	Add per expander
A	Waist	206	222	16
B	Hips	212	227	15
C	Upper Thigh	130	143	13
D	Knee	125	138	13
E	Ankle/Heel	104	114	10
F	End	87	95	8

*Provided for illustrative purposes only. Location of patient landmarks will vary depending on height and leg length. Contact your Lympha Press® dealer for assistance.

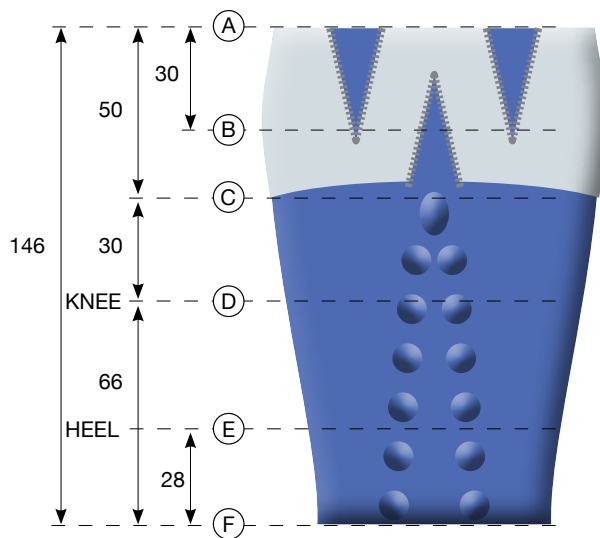

Use of the LymphaPod® is not recommended in the presence of one or more of the following conditions:

Known or suspected deep vein thrombosis or pulmonary embolism; during the inflammatory phlebitis process; acute infection of the affected limb, decompensated cardiac failure; severe arteriosclerosis or other ischemic vascular disease; any circumstance where increased venous and lymphatic return is undesirable. Due to movement of fluids in the body when using the device, be careful with patients who have heart disease. High pressure is not recommended for patients who have peripheral occlusion disease. The abdominal area should not be treated during pregnancy. If the LymphaPod® is to be used during pregnancy, the air cells over the abdomen must be disconnected prior to use.

Caution: Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner.

Distributed by: Lympha Press USA Ltd. 232 Park Avenue, Manalapan, NJ 07726

Tel: 1-888-596-7421, 732-792-9677, Fax: 732-792-9745

info@lymphapress.com www.lymphapress.com