



Lympha Press USA

# LYMPHA PANTS™



Note: We recommend Lympha Press® leg sleeves for treatment of uncomplicated unilateral lymphedema.

**Lower limb lymphedema can be associated with lymphedema in the torso, groin or abdomen.** Until now, manual massage has been the only option for torso therapy. Manual massage is time consuming and expensive. **For cases of severe lymphedema, manual massage is not enough.**

**Lympha Pants™ – torso treatment based on the principles of manual lymph drainage.** Designed in conjunction with leading manual lymph drainage therapists, the **lympa pants™** treat the entire lower body, including the groin and abdomen up to the chest.

### **Use the Lympha Pants™ for:**


- Lymphedema with groin, abdomen or torso swelling
- History of genital lymphedema
- Bilateral leg lymphedema
- Unilateral leg lymphedema with spillover to the healthy limb
- Intractable lower limb lymphedema that has not responded to other forms of treatment
- Elephantiasis

# LYMPHA PANTS™



## Features:

- 24 fully overlapping chambers treat both legs, groin, abdomen and lower back for maximal stimulation of all available lymphatics.
- Treats genital lymphedema, and prevents its development in patients who are susceptible.
- Bilateral treatment drains spillover, further enhancing lymph transport.
- Proven effective in hospitals, clinics, and homes worldwide.

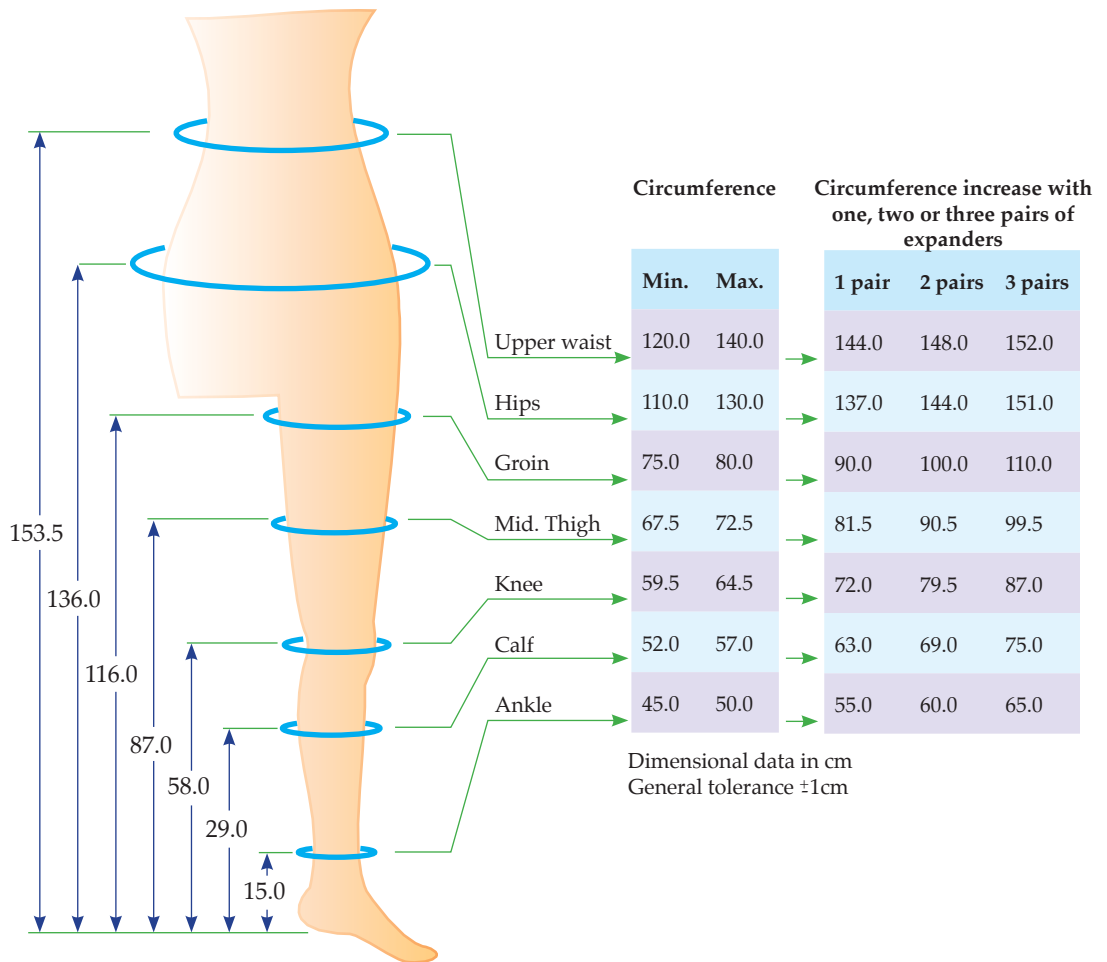

**The Lympha Pants™** accommodate wide variations in size. Our unique "custom boot" feature forms the sleeve to the patient's leg length. Expanders are available to increase circumference. Lympha Pants™ are also available with catheter access.

## What pump should I use?

**Lympha Press Plus®** - programmable lymph drainage massage

Lympha Press Plus offers five treatment programs, including our unique Pretherapy cycle based on the principles of manual lymph drainage.

**Lympha Press Mini®** - economical and effective

The Lympha Press Mini applies calibrated gradient pressure, ensuring that pressure at the torso is lower than on the legs, for comfortable and safe treatment. The Mini is ideal for home use, and is covered by most insurance plans.

