




Lympha Press USA

LYMPHA JACKET™

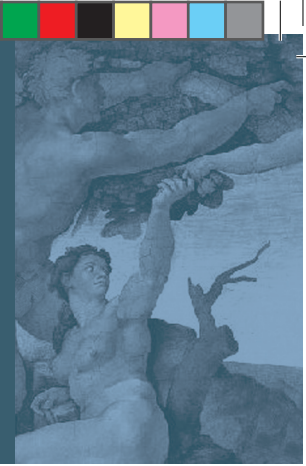

Complete treatment for upper torso lymphedema



Sequential pneumatic compression therapy with Lympha Press® is effective treatment for arm lymphedema. However, upper extremity lymphedema can appear in the torso as well as the arm, especially at the shoulder, upper back or chest.




Based on the principles of manual lymph drainage, the Lympha Jacket™ has 24 overlapping chambers, which compress and release in sequence to stimulate maximum lymphatic uptake. The Lympha Jacket™ treats the entire upper body, including both arms, shoulders, chest, back and abdomen. Bilateral treatment further enhances drainage and prevents spillover to the unaffected limb.

Use the Lympha Jacket™ for:

- Truncal lymphedema, including edema of the abdomen, chest wall, shoulder or back.
- Upper extremity lymphedema with shoulder, back or chest edema.
- Bilateral arm lymphedema.
- Unilateral arm lymphedema, if there is evidence of spillover to the unaffected limb.

Note: We recommend Lympha Press® arm sleeves, with or without shoulder attachment, for treatment of uncomplicated unilateral arm lymphedema.

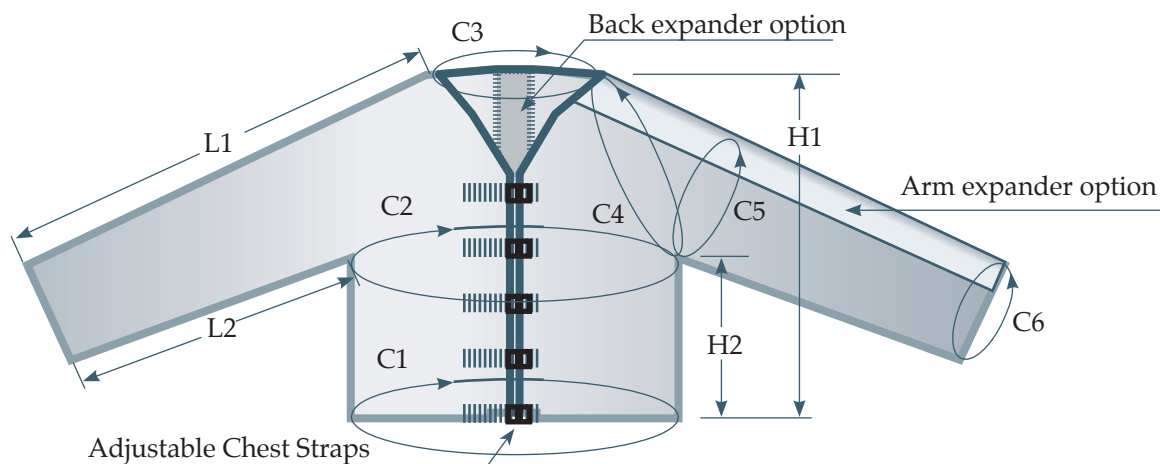


Jacket	Part No.	L1	L2	H1	H2	With chest straps				With back / arm expanders									
						C1&C2		C3		C4	C5	C6	C1&C2		C3		C4	C5	C6
						min	max	min	max				min	max	min	max			
Short arm	L4MJ2125	80	57	60	29	105	123	48	64	68	48	40	115	133	58	74	77	58	50
Long arm	L4MJ2135	96	72			105	123	48	64	68	48	40	115	133	58	74	77	58	50

All dimensions are in cm (centimeters) / General tolerance: ±1cm



Features:

- Complete upper body treatment, including both arms, shoulders, back, chest and abdomen.
- 24 fully overlapping chambers for superior directionality of flow, without gaps.
- Bilateral treatment prevents spillover and further enhances lymphatic transport.
- Fully adjustable across back, chest and arms. Expanders are available for an increase in size.

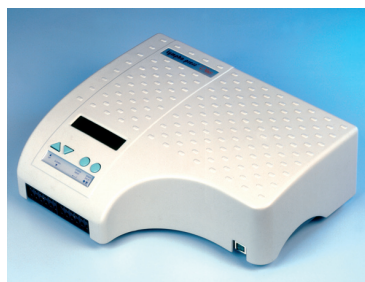

The Lympha Jacket™ accommodates wide variations in size. Our integral expander creates an exact fit across the chest. Additional expanders are available for an increase in size at the arms and back.

What pump should I use?

Lympha Press Plus® - programmable lymph drainage massage
 Lympha Press Plus offers five treatment programs, including our unique Pretherapy cycle based on the principles of manual lymphatic drainage.



Lympha Press Mini® - economical and effective
 The Lympha Press Mini applies calibrated gradient pressure, ensuring that pressure at the torso is lower than on the legs, for comfortable and safe treatment. The Mini is ideal for home use, and is covered by most insurance plans.

Lympha Press USA Ltd. 232 Park Avenue, Manalapan, NJ 07726
 Tel: 1-888-596-7421, 732-792-9677, Fax: 732-792-9745
 info@lymphapress.com www.lymphapress.com

